VELFORM®
Contour Sonic

Instruction Manual

Manual de instrucciones

Mode d'emploi

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





Manuaalinen

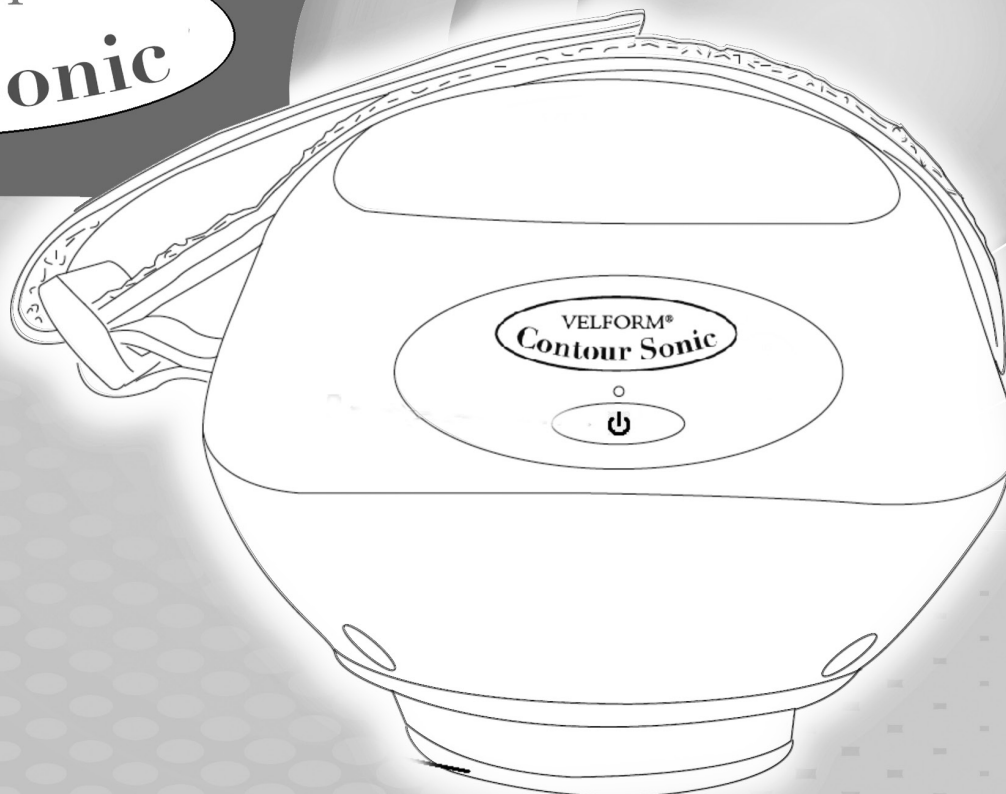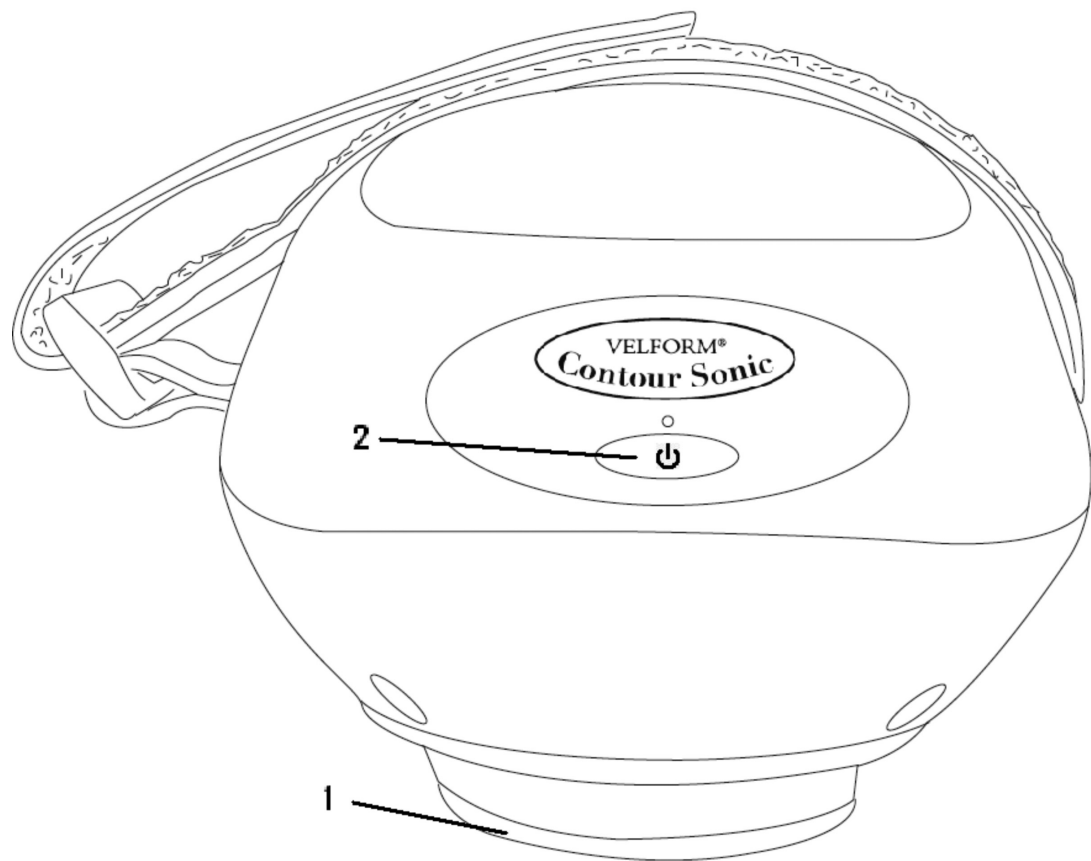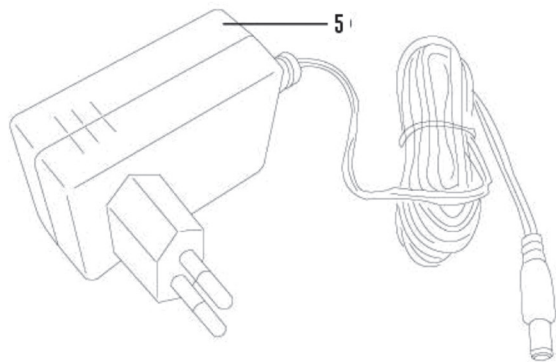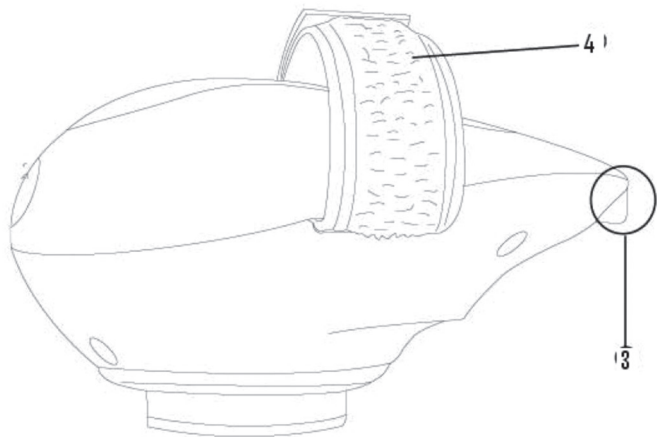A

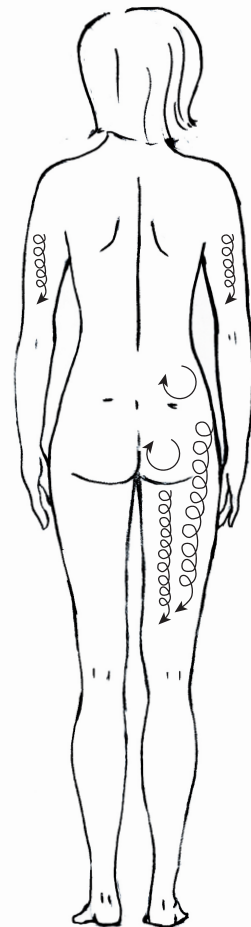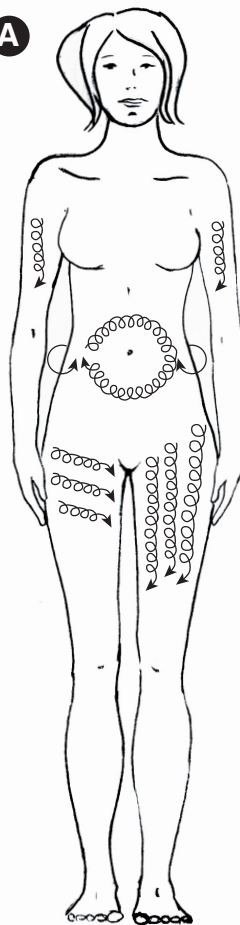VELFORM® CONTOUR SONIC INSTRUCTION MANUAL

Product description:

1. Head
2. Ultrasound button (ON/OFF)
3. Adaptor plug
4. Adjustable strap
5. Adaptor

SAFETY ADVICE AND WARNINGS

- Please read this instruction manual before using your Velform® Contour Sonic and keep it for future reference.
- Non-compliance with these instructions may lead to accidents.
- Do not stop moving the appliance and do not hold it in one position during application.
- Do not use the appliance if you are or think you may be pregnant.
- The appliance must not be used by people with heart conditions.
- The appliance must not be used by people with skin conditions, rashes, burns, scars, hyperpigmentation or broken capillaries.
- The appliance must not be used by people with silicone implants.
- The appliance must not be used by people who have undergone surgery in the preceding six months.
- When using Velform® Contour Sonic, we recommend applying plenty of Velform® Contour Sonic conductive gel on the area to be treated.
- Do not use the appliance on the face, neck, scalp, throat, shoulders, knees, genital area or spine. Avoid contact with the joints.
- If any type of allergic reaction occurs after using the product, stop using it immediately and consult a medical professional.
- Do not use the product on the same area for any longer than 10 minutes every two days.

IMPORTANT: Do not use the product if you have implants in the area to be treated, such as medical implants, an IUD, a prosthesis, screws, etc.

Use and care

- Clean the head with a handkerchief or cloth after each use.
- This appliance has been designed for domestic use only. Under no circumstances should it be used for commercial or industrial purposes.
- This appliance is not a toy. Children should be supervised to ensure that they do not play with it.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- This appliance is designed to be used by adults. It should not be used by anyone unfamiliar with this type of product, people with disabilities or children.
- Keep out of the reach of children and/or people with disabilities.
- Do not expose the appliance to extreme temperatures or strong magnetic materials.
- Do not use the appliance on pets or animals.

Instructions for use:

Before using Velform® Contour Sonic for the first time, remove the protective film from the head (1).

Connect the adaptor (5) to the adaptor plug (3) and connect the adaptor to the power supply.

Place your right or left hand, depending on which is more comfortable, on the back of the appliance and use the adjustable strap (4) to ensure it is securely fastened.

When using Velform® Contour Sonic, we recommend applying Velform® Contour Sonic conductive gel on the area to be treated.

To turn on the device, press the ON/OFF button (2). If you have done it correctly, the green pilot light will light up.

Next, place the head on the area to be treated and move it as shown in Figure A. Do not stop moving the appliance and do not hold it in one position during treatment.

Press the ON/OFF button again to turn the device off.

The device is programmed to turn itself off after 10 minutes of operation. When the ON/OFF pilot light goes out, you will know that the session has ended. The appliance should not be used in the same area for more than 10 minutes. Do not repeat treatment in the same area the following day. You must allow the area to rest for a full day before treating it again.

Technical Specifications of Velform® Contour Sonic

Brand: Velform® Contour Sonic

Input: 12V  0.5A 6W

Frequency: 1 MHz

Polarity: 

Technical specifications of the adapter:

Model: GS2T-0120500-1

Input: 100-240V~50/60Hz 0.5A

Output: 12V  0.5A 6W

Polarity: 

 This appliance has a double insulating system, so it does not need ground connection. Please make sure that the tension of the mains corresponds to the value indicated on the appliance's information badge.

Warranty

This product is covered by a warranty against manufacturing defects subject to the time periods stipulated by the legislation in force in each country.

This warranty does not cover damages resulting from inadequate use, negligent commercial use, abnormal wear and tear, accidents or improper handling.

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CERTIFICATES AND DIRECTIVES

This product complies with Electromagnetic Compatibility Directive 2004/108/EC, RoHS Directive 2002/95/EC and Low Voltage Directive 2006/95/EC.

This product must not be disposed together with the domestic waste. This product has to be disposed at an authorized place for recycling of electrical and electric appliances. By collecting and recycling waste, you help to save natural resources, and make sure the products is disposed in an environmental friendly and healthy way.

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ISL GmbH
Best Direct International PTY LTD
Best Direct International Ltd,