

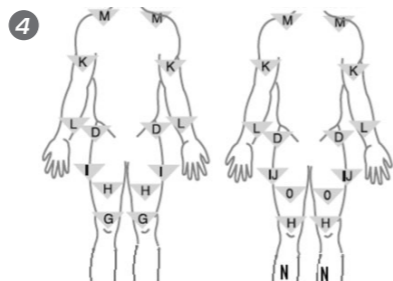
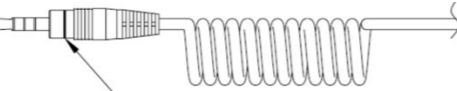
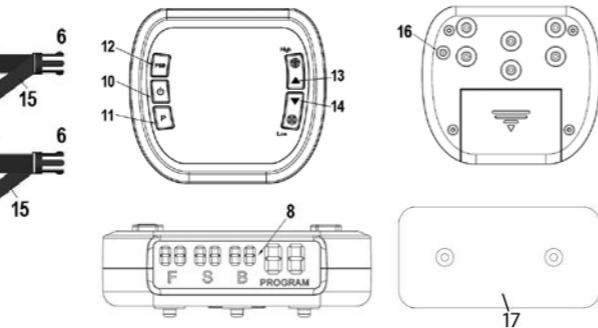
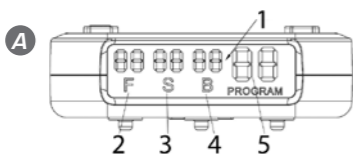
-  *Instruction Manual*
-  *Manual de instrucciones*
-  *Mode d'emploi*
-  *Bedienungsanleitung*
-  *Istruzioni per l'uso*
-  *Manual de instruções*
-  *Gebruiksaanwijzing*

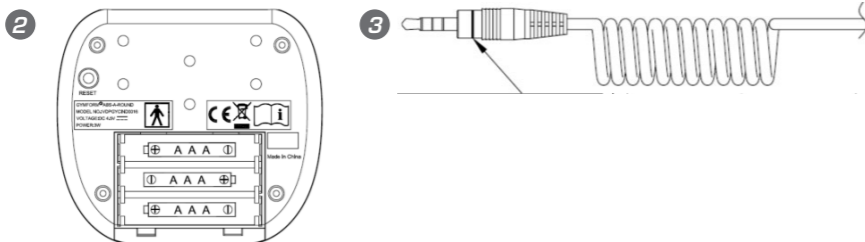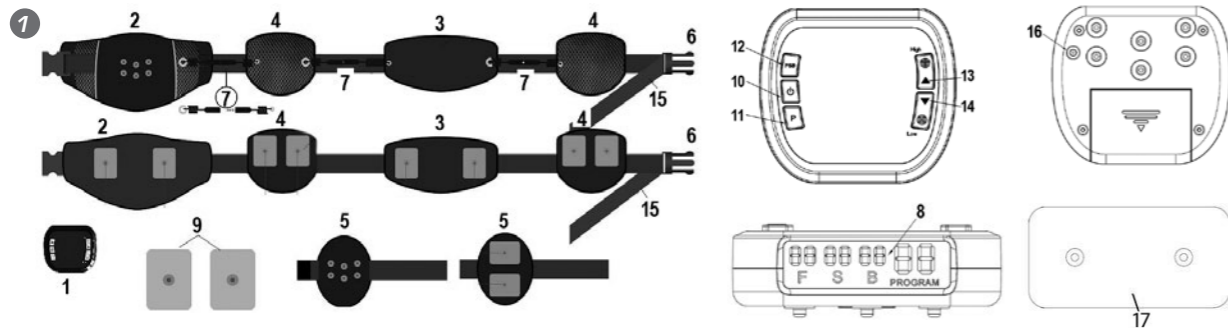Congratulations on your purchase of Gymform® ABS-A-ROUND. Your safety is our primary concern, so please ensure that you familiarize yourself with all the information in this manual before using your Gymform® ABS-A-ROUND. It is essential that anyone who intends to use this device reads and understands the instructions provided. An informed user is a safe user. Thank you for purchasing the Gymform® ABS-A-ROUND. We hope you enjoy using it.

WARNING

Before embarking on any exercise programme, including the Gymform® ABS-A-ROUND, consult your general practitioner and your physiotherapist. Ensure that you read the instructions, warnings and precautions carefully before using your Gymform® ABS-A-ROUND. Gymform® ABS-A-ROUND is designed for healthy adults. However, the device should be used with care. Read the following guidelines and the instruction manual before using the product.

- Applying electrodes in the thorax area may increase the risk of cardiac fibrillation.
- Using the Gymform® ABS-A-ROUND while connected to high-frequency surgical equipment may cause burns at the contact points of the Gymform® ABS-A-ROUND and may result in damage to your Gymform® ABS-A-ROUND.
- Operation of a short-wave therapeutic device or microwave oven in the general vicinity (e.g. 1 metre away) may cause instability in the output power of the Gymform® ABS-A-ROUND.
- This is a powerful electrical muscle stimulation device and may cause injury if you start with the most advanced exercise programme right away. Start with a lower level and gradually increase the intensity.
- Do not use on children or animals.
- Do not use while driving or sleeping.
- Do not use in the bathroom or in other places with a high amount of moisture or humidity.
- The batteries should be removed from the device before it is disposed of.
- The batteries should be safely disposed of.
- This device may be used by children aged 8 and over, by people with reduced physical, sensory or mental abilities, and by people lacking in experience and knowledge if they are appropriately supervised or given suitable instructions on the safe use of the device and understand the dangers it use entails. Children should not be allowed to play with the device. The cleaning and maintenance to be performed by the user should not be done by children without supervision.

DO NOT use Gymform® ABS-A-ROUND:

- If you have an electronic implant (such as a pacemaker) or suffer from another heart problem.
- If you are, or think you might be, pregnant.
- If you have cancer, epilepsy or are under medical supervision for a cognitive disorder.
- If you are breastfeeding.
- If you have a high fever.

Women: before using your Gymform® ABS-A-ROUND, you should wait:

- At least six months after the birth of your child (after consulting your doctor).
- A month after having an intrauterine device (IUD) fitted.
- At least three months after having a caesarean section (after consulting your doctor).
- Until after the first days of your menstrual period.

Before using Gymform® ABS-A-ROUND you should obtain permission from your doctor or physiotherapist:

- If you suffer from back problems. In this case, the device should only be used at low intensity levels.
- If you suffer from a serious disease or injury not mentioned in this section.
- If you were operated on recently.
- If you are diabetic and take insulin.

When putting on the belt, always remember:

- To place the large belt ONLY on your abdomen, obliques and lower back.
- To place the small belt ONLY on the areas indicated below.
- Not to place the device on your neck, on your chest or on your back in such a way that electricity passes through your heart (e.g. an electrode on your chest and another on your back), genitals, scars moles, varicose veins or head).
- Not to place the belt on recent scars, inflamed skin, infected areas or areas that are susceptible to acne, thrombosis or other vascular problems (e.g. varicose veins), or any part of your body with limited sensitivity.
- Not to place the belt on areas with injuries or limited mobility (e.g. fractures or sprains).
- Not to place the belt directly on metal implants. It should be placed on the closest muscle.

Possible Adverse Reactions

There have been isolated reports of skin reactions caused by the use of electro-stimulation devices, including allergies, prolonged skin redness and acne.

If any of these symptoms occur, switch off the device and consult your doctor.

If you have never exercised before or have not exercised for a long time:

Consult your general practitioner.

Precautions:

- Do not use Gymform® ABS-A-ROUND while you are drowsy or asleep.
- Do not use Gymform® ABS-A-ROUND if your body is in contact with water.
- Keep the device out of reach of children.

- Do not use Gymform® ABS-A-ROUND in combination with any other device that transmits electrical current to your body or with any other machine.
- Do not connect the metal snap fasteners or electrodes to any other object.
- Do not modify the device.
- Use only the electrodes supplied by Industext.
- Do not attach two electrodes together.
- If any of the electrodes is not attached properly, you will notice a itching sensation.

Note:

- The product may malfunction as a result of electrostatic disturbances. If this occurs, simply press the RESET button (16) to resume normal operation.

CONTENTS:

- 1 removable Gymform® ABS-A-ROUND unit with six programmes
- 1 Gymform® ABS-A-ROUND belt
- 1 mini belt for working the leg and arm muscles
- 6 self-adhesive electrodes (2 of which are double electrodes)
- 1 instruction manual
- 1 sports bag
- 1 diet booklet

PRODUCT DESCRIPTION (Fig.1):

1. Gymform® ABS-A-Round unit
2. Belt: abdominal part
3. Belt: lower-back part
4. Belt: abdominal-oblique part
5. Mini belt: arm and leg part
6. Side-release buckle
7. Belt connectors
8. LCD display
9. Electrodes for abdominal and lower-back areas
10. ON/OFF button
11. Programme-selection button
12. Button to select belt part
13. Button to increase intensity
14. Button to reduce intensity
15. Belt adjustment strap
16. RESET button
17. Double electrodes for abdominal-oblique areas

Putting on the Belt:

Set-up:

Open the battery compartment cover (batteries not included) on the back of the unit (1) and insert three AAA batteries. Insert the three batteries in accordance with the polarity marks. The polarity is indicated on the battery compartment for proper placement, as shown in the illustration. Replace the battery compartment cover. (Fig. 2)

Putting on the Belt:

Attach the Gymform® ABS-A-ROUND unit (1) to the front of the belt using the six snap fasteners. Ensure that they are all well fastened or the belt will not work properly.

Attach two electrodes (9) to the abdominal part (2) and two electrodes to the lower-back part (3) of the belt.

Attach two double electrodes (17) to the two parts of the abdominal and obliques (17).

Make sure the belt connectors (7) between the abdominal, lower-back and oblique parts of the belt are properly connected by ensuring that the indicated mark on each belt connector (7) is not visible. (Fig. 3)

Remove the protective plastic backing from the electrodes (9) on the abdominal part of the belt (2). Place the electrodes on your abdominal area and press to stick down. Then remove the protective plastic backing from the double electrode (17) on the left oblique part of the belt and place the electrodes on your oblique muscle, which is located above the hip bone. When the left oblique part is in position, adjust the belt so that it is taut between the abdominal (2) and oblique (4) parts. Repeat this procedure to attach the electrodes of the lower-back part (3), followed by the right oblique part. Ensure that the belt is properly adjusted between all parts so that the electrodes are in contact with the skin. Then fasten the belt using the side-release buckle (6).

When you have finished using the belt, to ensure the electrodes do not stick to the skin, remove them one at a time by pulling them from the edge and then replace the protective plastic backing on the electrodes. This will help extend the useful life of the electrodes.

Using the Belt

Press the ON/OFF button for one second to switch on your Gymform® ABS-A-ROUND.

Figure A:

1. LCD SCREEN
2. Abdominal part
3. Obliques part
4. Lower-back part
5. Program

Programme Selection

Programme 1	General massage
Programme 2	Fine tuning

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Best Direct International Pty LTD
Best Direct International Ltd,